

## Installation Guide


### Philips FMS/MMS Module Lock

The purpose of this guide is to describe installation of the Module Lock.

#### Parts Reference

The following parts and hardware are included in this installation kit (hardware not shown):

Item #	Description	Qty
1	Main Enclosure	1
2	Support Bracket	1
3	#8-32 x 1/4" Button Head Machine Screw, Tamper-proof	4
4	5/32" Hex Wrench, Tamper-proof	1

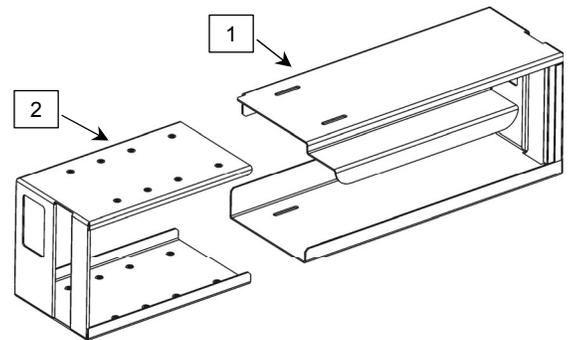

**Tools required:** 5/32" tamper-proof hex wrench (provided).

**Installation Note:** The Support Bracket is designed to accommodate an FMS (handles on both ends) or with MMS Extension (handle on one end). Installation is shown only on FMS with MMS Extension.

#### Installing the Module Lock

1. With required modules installed, slide FMS into Main Enclosure as shown below. **Installation Note:** It is not necessary to remove a previously-installed mounting clamp from the FMS to install the Module Lock.

