

# PHILIPS

IPL – Hair removal device

Lumea Advanced

For body and face

With skin tone sensor

15 minutes to treat lower legs

+ Satin Compact pen trimmer



BRI921/00



## Be hair free for at least 3 months\*

### At least 85% hair reduction in 3 treatments\*\*

Philips Lumea Advanced works wonders to prevent the reappearance of hair on the body. Gentle pulses of light, applied regularly, keep skin silky smooth every day.

#### Long-lasting hair removal with IPL technology

- Expert IPL technology at home, developed with dermatologists
- Suitable for a wide variety of hair and skin types
- FREE Lumea App for a personal coach at your fingertips

#### Safe, fast and effective on your body and face


- Proven safe and effective treatment
- Integrated skin tone sensor

#### Specialised treatment of different body areas

- Safe and effective even on sensitive areas
- Large treatment window for fast application
- Precision attachment for extra safety on the face
- Special Beauty Edition including Satin Compact pen trimmer


# Highlights

## Long lasting smoothness

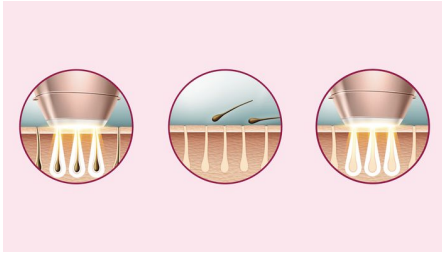

IPL stands for Intense Pulsed Light. Philips Lumea applies gentle pulses of light to the hair root, putting the follicle into a resting phase. As a consequence, the amount of hair your body grows gradually decreases. Repeating the treatment leaves your skin beautifully hair free and touchably smooth. Treatment to prevent hair regrowth is safe and gentle, even on sensitive areas. Philips Lumea is clinically tested and developed with dermatologists for easy and effective treatment, from the comfort of your own home.

## Proven effectiveness

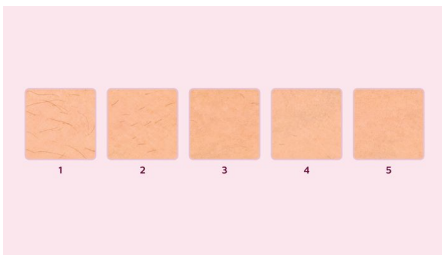


Objective studies show at least 85% hair reduction in as little as 3 treatments\*\*. The first 4 treatments must be carried out every 2 weeks, after which you should already be able to see the desired results. To maintain your results, simply touch up every 4 weeks. After just 8 touch-up treatments, you can be hair-free for at least 3 months\*. To enlarge this picture, please click on the image in the image gallery at the top of this page.

## Suitable hair and skin types



Philips Lumea Advanced works effectively, easily and safely on a wide variety of hair and skin types. It works on naturally dark blonde, brown and black hair and on skin tones from very white to dark brown. As with other IPL-based treatments, Philips Lumea cannot be used to treat white/grey, light blonde or red hair and is not suitable for very dark skin. This is due to the high contrast required between the pigment in the hair colour and the pigment in the skin tone.

## Skin Tone Sensor

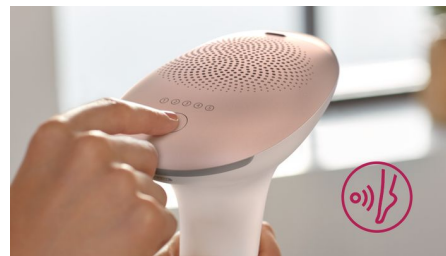

The integrated skin tone sensor measures the treated skin complexion at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for this product, it will automatically stop emitting pulses.

## Application areas



Use with confidence to remove hair on face (upper lip, chin and sideburns) and body parts including legs, underarm, bikini area, stomach and arms.

## Fast treatment on legs




Big treatment window for fast application on large areas like legs.

## Extra safety for the face



With the additional integrated light filter, it can be used for safe treatment of facial hair on the upper lip, chin, cheeks and other sensitive areas.

## FREE Lumea App

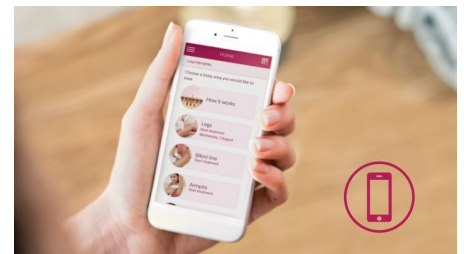

The Lumea App is free to download and offers a unique and personalised IPL experience – just for you. It is your very own personal coach to ensure you are getting the most out of your Lumea and using it in the correct way to achieve long-lasting results. The App helps you create your personalised treatment schedule per body area, with tips and advice during each treatment. All you need to do is check the app regularly to make sure you are up to date with any notifications or reminders.

# Specifications

## Application areas

**Body areas:** Bikini, Underarms, Arms, Legs, Belly

**Face areas:** Chin, Sideburns, Upper lip

## Attachments

**Body attachment (4 cm<sup>2</sup>):** For use below the neckline

**Face attachment (2 cm<sup>2</sup>):** For use below the eyes

## Application mode

**Slide and Flash:** For easy motion application

**Corded / cordless use:** Corded

**Stamp and Flash:** For treatment on curvy areas

## Application time

**Bikini:** 4 min

**Underarm:** 2.5 min

**Face areas:** 2 min

**Lower leg:** 15 min

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**Skin Tone sensor:** Detects your skin tone

## Power

**Voltage:** 100–240 V

**Corded use**

## Items included

**Instructions for use:** User manual

**Storage:** Pouch

**Adapter:** 24 V/1500 mA

**Accessory:** Cleaning cloth

## Service

**Warranty:** 2 years global warranty

## Pen trimmer

**Ease of use:** Compact handle

**Accessories:** Eyebrow comb, Cleaning brush

**Power:** AAA battery

